
MEDI GUN

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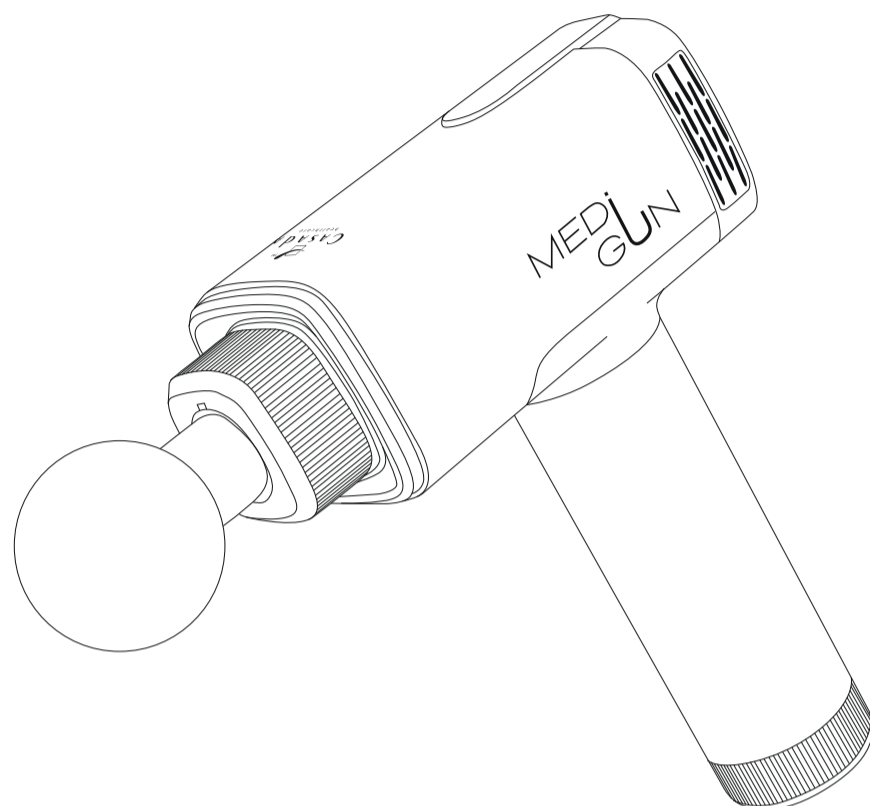Casada International GmbH
www.casada.com

DAS PRODUKT

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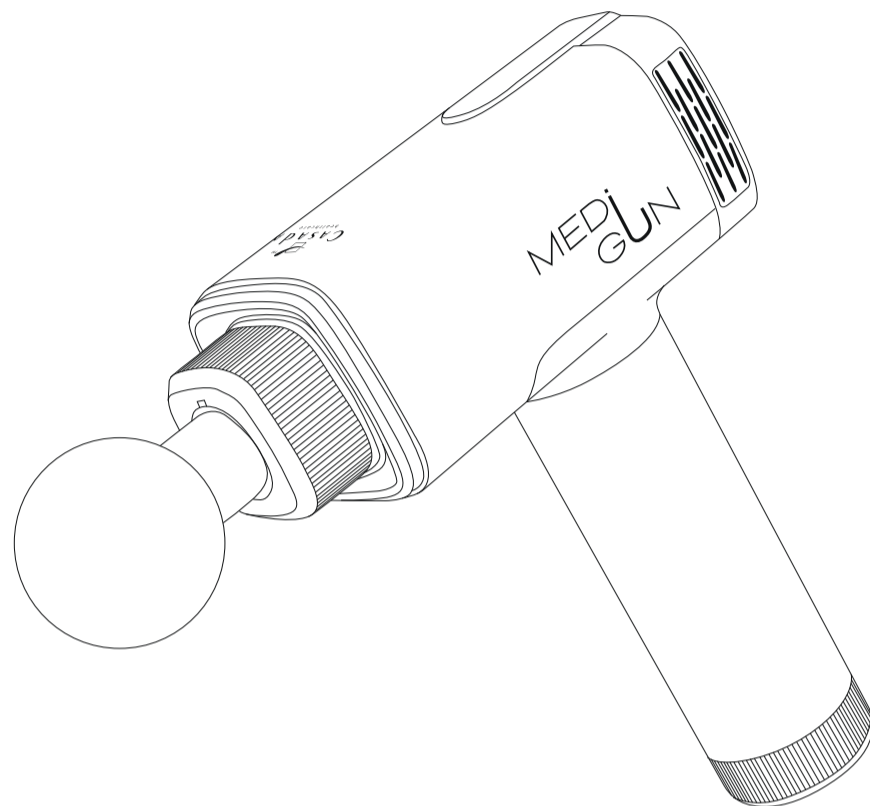
Casada International GmbH
www.casada.com

PRODUCT INTRODUCTION

MediGun by Casada is a wireless, state-of-the-art massage device that helps to reduce muscle soreness and tension, promote blood circulation, optimize mobility and improve the general condition of the body's musculoskeletal system.

CORDLESS / RECHARGEABLE
24V Lithium-Ion Battery

WEIGHT
2.2lbs / 1.0kg

SPECIFICATIONS | TECHNICAL DATA

SPECIFICATIONS

- 1 Advanced technology with powerful brushless High-Torque Motor
- 2 Interchangeable head attachments (x4) – Fork, Bullet, Flat, Round
- 3 Rechargeable Lithium-Ion battery (More than 2 Hours using time when full charged)
- 4 Gear indicator
- 5 Speed Settings for Optimal Pressure
- 6 Battery indicator
- 7 AC interface
- 8 ON/OFF Power Switch

TECHNICAL DATA

Frequency: Level 1 – 30Hz / Level 2 – 40Hz / Level 3 – 53Hz
Weight: 2.2lbs / 1.0kg
Charge: 110–240V, 50/60Hz – 5A, 26V – DC – 1A

USAGE INSTRUCTIONS

⚠️ USAGE INSTRUCTIONS PLEASE READ PRIOR TO USE

CHARGING

- Please fully charge up to six hours when first use.
- When charging, connect the DC end of the supplied 26V adapter to the battery's charging port and plug AC side cable into wall outlet, then turn the power switch to the ON position to charge.
- The LED light flashes to indicate the battery level and indicating it is charging.
- When all LED stop flashing and remain lit, it indicates that the battery was fully charged.
- Battery can be recharged at any time and at any battery level.
- It is not recommended to fully deplete the battery to the red LED level.
- This device average run-time is about 2 hours depending on speed level and pressure applied during using.

OPERATING

1. Turn the power switch on the bottom of the battery to the "ON" position to start the power. Turn the power switch to the "OFF" position to turn off the power.
2. With the power on, short-press the gear button once to turn device on to level one, twice to turn device to level two and three times to level three. Pressing the gear button again will turn the device off. The blue LED lights would indicate the corresponding frequency level.
3. Massage desired body part at the desired speed applying pressure as needed.
4. With power off, you can change the massage head and pull it out by hand.

MAINTENANCE, CLEANING & STORAGE

- Use a slightly damp towel to clean device and a dry with a soft cloth.
- For storage or travel, remove battery and store in supplied carton box.

DISCLAIMER

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 30 minutes before re-using.
- Unplug the charger from device battery after charging or prior to use.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

SAFETY INSTRUCTIONS

DANGER
To reduce risks of electric shock, fire, and personal injury, this product must be used in accordance with the following instructions:

- For adult use only.
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the MediGun only on the body's soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
- Use only the applicators that provide the best desired effect.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the MediGun.
- Do not immerse in water or allow water to enter the ventilation holes of the MediGun.
- Do not drop or misuse the MediGun.
- Only recharge with supplied 24VDC MediGun charger.
- Carefully examine device and battery before each use.
- Do not tamper or alter MediGun in any way.
- Never leave the MediGun operating or charging unattended.

DISCLAIMER

Please do not use the MediGun, or any vibration device without first obtaining approval from your doctor if any of the following apply:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylolysis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health.

Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first. Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

WARRANTY

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt. The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period. The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorised dealer. The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product. If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer. Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medical treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Congratulations!

With the purchase of this massage device you have shown awareness of your health.

Since being established in 2000, the name CASADA stands for high quality products on the healthspa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide. Constant further development, an unmistakable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.